

Reduced Viscosity Analysis of Polyurethane

Application Note

Materials Testing and Research, Polymers

Authors

Greg Saunders, Ben MacCreath
Agilent Technologies, Inc.

Introduction

Polyurethanes are soluble in several medium to high polarity solvents. In DMF, the temperature is increased to reduce eluent viscosity so reducing operating pressure. Agilent PLgel 5 μm MIXED-C columns are well suited to the analysis of polyurethanes.



Agilent Technologies

PLgel 5 μ m MIXED-c columns are designed for rapid polymer analysis. With its linear calibration up to 2 million MW, this is the column of choice for highest resolution and accuracy in molecular weight distribution analyses. Rapid solvent change capability, excellent temperature stability and the high resolution of the PLgel 5 μ m MIXED-C also provide the versatility essential for today's R&D laboratory.

Conditions

Columns:	2 x PLgel 5 μ m MIXED-C, 300 x 7.5 mm (part number PL1110-6500)
Eluent:	DMF + 0.1% LiBr
Flow Rate:	1.0 mL/min
Temperature:	60 °C
Detection:	390-MDS Multi Detector Suite (differential refractive index)

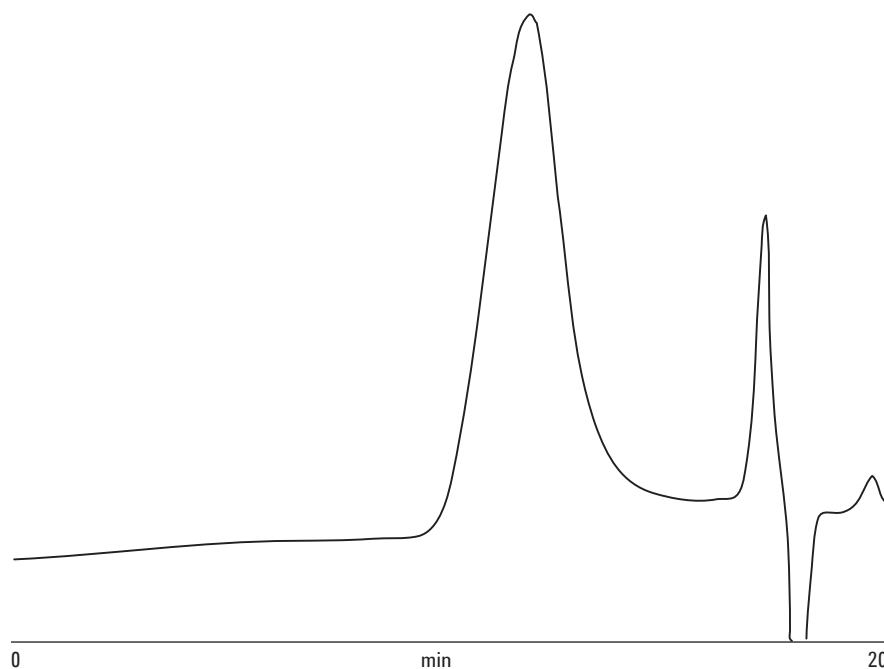


Figure 1. Analysis of polyurethane using PLgel 5 μ m MIXED-C columns

www.agilent.com/chem

This information is subject to change without notice.

© Agilent Technologies, Inc. 2010

Published in UK, September 29, 2010

SI-1823



Agilent Technologies